# USER'S MANUAL

---

MNPG314 Rev.0 del 22/05/2020

Infrared Forehead Thermometer

# FI06

**ATTENTION: for correct measurement of the temperature, the thermometer and the patient must have been at the same room temperature for 30 minutes. Sweat on the forehead or just wearing a cap can lead to incorrect measurements!**

Manufacturer

**Manufactured by: Ningbo Shangcun Electronics Co.,LTD**  
No.9 Yongning Rd, Yuyao, Zhejiang, China  
Tel. (86) 755 2765 2471 • Fax (86) 755 2765 2674

**GIMLE Kamila Sienkiewicz**  
**Ul. Sw. Judy Tadeusza 108, 52-116 Iwiny, Poland**  
Tel. +0048 601 669717 • Fax. +0048 601 669717

**Distributed by: I.A.C.E.R. S.r.l.**

Via S. Pertini, 24/a • 30030 Martellago (VE)  
Tel. 041.5401356 • Fax 041.5402684

Ningbo Shangcun Electronics Co., LTD is a Chinese manufacturer of medical devices (certified CE n°0197 from the Notified Body n°0197 TÜV Rheinland LGA Products GmbH).

### Product introduction and classification

The FI06 device has the following classification:

- class IIa (Directive 93/42/CEE, Annex IX, rule 10 and further amendments);
- Type B applied part.

The infrared thermometer is a thermometer that measures the temperature of human body,

quickly and accurately, by using the principle of receiving infrared.

### Product features

- Specially designed for measuring the body temperature.
- Has temperature alarm sound indicators when this exceeds the default threshold value.
- Capable of storing 32 sets of measurement data.
- Colours backlight LCD digital display (red, orange, green).
- Measurement in Celsius/Fahrenheit.
- Providing the function of auto shutdown to save electric power.
- Selecting the body temperature measurement mode to measure the temperature of any object with its temperature less than 100°C and its emissivity equal to 0.95.

### Technical parameters

Characteristic	Specific
Power supply	DC 3V 2AAA alkaline Battery
Power consumption	When off $\leq 10\mu\text{W}$ when measurement $\leq 30\text{mW}$
Dimensions (length x width x height)	125 x 74 x 38mm


<b>Characteristic</b>	<b>Specific</b>
Weight	82gr (without batteries)
Body measurement range	32.0°C-42.9°C (85.9°F-109.2°F)
Body Accuracy	±0.2°C (±0.4°F)
Object measurement range	0°C-100°C (32°F-212°F)
Object Accuracy	±1.0°C (±1.8°F)
Resolution	0.1°C/°F
Power level indicator	Indication for low power level (≤2.5W)
Memory	Automatic recall of 32 memories
Display	LCD back-lit display
Reading scale	Celsius o Fahrenheit
Automatic shut off	in 20 seconds
Standards	EN60601-1 ISO80601-2-56 ARTS 1965-98
Batteries	2AAA alkaline Battery

Characteristic	Specific
Working condition	Temperature: 16-35°C
	Humidity: ≤85%
Storage condition	Temperature: -20-55°C
	Humidity: ≤93%

### Contents of the packaging

- N°1 Thermometer;
- N°2 AAA batteries;
- N°1 user's manual.
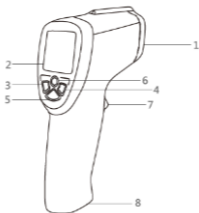
### Device description and comands

- (1) IR Sensor
- (2) LCD display
- (3) Memory inquiry or Minus
- (4) Memory inquiry or Add


- (5) C/F Chang or Advanced setting
- (6) ON/OFF button
- (7) Measurement trigger
- (8) Battery cover

**Function definition of all icons**

Function	Icon	Details	
Battery level		When it is visible	The battery is in low level, but the thermometer is functional properly. Please replace battery asap.
		When it flashes	The battery is exhausted, and thermometer cannot function properly. Please replace battery immediately.
		When it is invisible	Battery power is sufficient
Laser		When it is visible	Laser on (when the machine installed laser head)
		When it is invisible	Laser off (when the machine installed laser head)
Measurement mode	<b>Body temp</b>	Body mode	
	<b>Object temp</b>	Object mode	

Function	Icon	Details
Reading scale	°C	Celsius reading
	°F	Fahrenheit reading
Reading display	<b>8888</b>	Temperature value
Memory	<b>M</b>	Memory value
Buzzer		Buzzer on
		Buzzer off
Bluetooth		Bluetooth on/off (optional)

## Function definition of all buttons

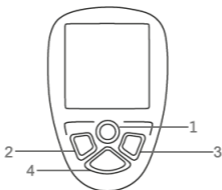N.	Description
1	Press shortly: Backlight ON/OFF Press for 3 seconds: Object/Body temperature change
2	Memory inquiry
3	Memory inquiry

N.	Description
4	Multifunction button Press shortly: voice ON/OFF Press for 3 seconds: function setting

## Setting

Menu	Function	-	+	Default	Remarks
F1	Deviation value setting	Decrease 0,1°C	Increase 0,1°C	0.0°C	Object mode is useless
F2	Alarm value setting	Decrease 0,1°C	Increase 0,1°C	38°C	Effective range: ±3
F3	°C/°F	°C	°F	°C	

### NOTICE

1. Temperature under human body mode is obtained from dynamic compensation of environmental temp and forehead surface temp.
2. Object temperature mode is to test surface temperature of an object. The temperature get from forehead under this mode is merely temperature of forehead surface but not body temperature.
3. Deviation setting can adjust measurement value from -3.0°C to 3.0°C according to testing distance, surrounding temperature, skin difference, etc. Defaulted value is 0.0°C.

**Measurement****BODY TEMPERATURE**

- Press measurement button to turn on thermometer and it displays boot screen (fig. 9.1). After POST and two beeps, it will display value of room temperature and be ready for measurement (fig. 9.2).
- Make sure the thermometer is under body mode.
- Keep distance at less than 5cm from upper eyebrows to the probe (fig. 9.3). Press measurement button and when it gives a “beep” measurement value is exceeding alarm value (Defaulted value is 38°C), it will give 3 beeps as an indication.
- After measurement, if the thermometer is idle in 20 seconds, it will display “OFF” (fig. 9.5) and gives a beep and shut off automatically.

(fig9.1)

(fig9.2)

(fig9.4)

(fig9.5)

(fig9.3)

**NOTICE**

1. Keep the sensor and probe cavity clean before and after use.
2. To ensure the accuracy of measurement, it is recommended to start measurement after ten minutes when carrying the thermometer to a new environment.
3. Wait for 10 minutes to measure body temperature after measuring extremely high temperature or extremely low temperature objects.
4. Wait for 5 minutes to start a measurement when a measuring target (object or human) is from an environment with enormous difference in temperature.
5. Breeze, water, sweating, cosmetic on forehead may affect measurement. Do not measure body temperature in 30 minutes after exercise, bath, or meals.

**OBJECT TEMPERATURE**

- Press measurement button to turn on the thermometer (fig.9.6).
- After measured room temperature make sure the thermometer is under object mode (fig.9.7-9.8).
- Keep vertical distance at 5cm to 15 cm from object to measurement probe. Press measurement button and when it gives a beep measurement is finished and value is idle in 20 seconds. It will display "OFF" (fig. 9.9) and gives a "beep" and shut off automatically.

(fig.9.6)

(fig.9.7)

(fig.9.8)

(fig.9.9)

The value under this mode is object surface temperature instead of core temperature.

### Exceeding measurement range

Body mode:

- When measurement value is lower than 32.0°C, it displays Lo (fig.9.10) and gives 4 beeps.
- When measurement value is higher than 42.9°C, it displays Hi (fig.9.11) and gives 4 short beeps.

Object mode:

- When measurement value is lower than 0°C, it displays Lo (fig.9.12) and gives 4 beeps.
- When measurement value is higher than 100.0°C, it displays Hi (fig.9.13) and gives 4 short beeps.

**ATTENTION:** When surrounding temperature is lower than 10.0°C or higher than 40.0°C, it displays Err. It is not allowed to measure, or accuracy is not assured.

(fig9.10)

(fig9.11)

(fig9.12)

(fig9.13)


(fig9.14)

### Battery replacement

- Open the battery lid and take out exhausted battery.

- Put into 2AAA alkaline batteries and close battery lid. After new battery is installed, thermometer will give 2 beeps. If there is no beeps, check if the positive and negative pole is correct (fig.9.1).

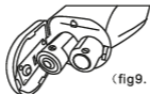(fig9. 1)

**ATTENTION:**

- Take out battery in case the thermometer is not used for long period. Do not put the battery to fire.
- Do not throw old or defective batteries into household waste, for disposal use the appropriate containers or otherwise according to current regulations (WEEE). Keep the batteries out of reach of children.

**Maintenance**

- Make sure the sensor and probe cavity are clean otherwise it will affect accuracy. Cleaning method for probe:
  1. Use the cotton stick or soft cloth with water or alcohol to wipe the casing.
  2. Use the cotton stick or soft cloth with alcohol to wipe the sensor surface or probe cavity gently. Do not use thermometer before alcohol is vaporized.

- Read this manual book thoroughly before use. Make sure battery is well installed.
- It is not allowed to put the thermometer in any liquid or expose to strong sunlight or extremely low temperature.
- Strong crash or hit to the product will cause its damage.
- Do not dismantle this thermometer by yourself.
- Keep the thermometer from children's reach.
- Do not use the thermometer under circumstance of strong electromagnetic interfere.
- The measurement results are probably fluctuating due to improper measurement ways. Please practice adequate measurements to improve your skill. **The measurement results can not supersede a doctor's diagnosis.**
- If the thermometer should not work contact the distributor or I.A.C.E.R. Srl.

#### **INFORMATION ON DISPOSAL**

To avoid contaminating or infecting personnel, the environment or other equipment, ensure that the equipment is properly disinfected or decontaminated before it is disposed of in accordance with the law of its country for equipment containing electrical and electronic components.

**Trouble shooting**

<b>Description</b>	<b>Solutions</b>
LCD display "Lo" or "Hi"	<ol style="list-style-type: none"><li>1. Breeze, water, sweating, cosmetic on forehead may affect measurement.</li><li>2. Check deviation value setting. Defaulted value is 0.8°C.</li><li>3. While if the testing environmental temp changes so enormously or if the thermometer is used directly from high-temp object to very low-temp one, the measurement difference will happen. The thermometer should be kept in a relative stable environment for 10 minutes to get heat balance before starting a new measurement.</li><li>4. Ensure measurement distance is 5cm to 8cm.</li></ol>
No response when pressing measurement button.	<ol style="list-style-type: none"><li>1. Take out and reassemble battery.</li><li>2. Check if the thermometer is under menu setup. In procedure of menu setting, thermometer is unable to measure and therefore no response.</li></ol>
No display or improper display	Take out battery and install battery again.
No beeper	Check if the beeper is switched off.

Description	Solutions
Shut off right after switching on	Check battery level or take out and install the battery again.

### Warranty, Support and Spare parts

FI06 is covered by 2-year warranty from the date of purchase on the electronic parties. The warranty does not cover parts subject to normal wear and tear of use (batteries) and all parts that are defective due to negligence or carelessness in use, incorrect maintenance or in the event of equipment tampering and intervention on the equipment by personnel not authorized by the manufacturer or authorized dealer.

If you have any questions about the device or this warranty, please contact I.A.C.E.R. Srl I TECH MEDICAL DIVISION. The distributor is the one and only allowed to operate with technical assistance. For any technical assistance contact the dealer or I.A.C.E.R. S.r.l.

Technical documentation related to repairable parts could be attached, but only with previous authorization from the company and only after giving proper training to the staff employed in technical assistance. The distributor makes available at any time the original spare parts for the equipment. Please contact I.A.C.E.R. S.r.l. In order to preserve the warranty, the functionality and the security and safety of the product, it is highly recommended to use exclusively the spare parts given by the manufacturer.

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- Si le thermomètre ne fonctionne pas, contactez I.A.C.E.R. Srl.
**Distributed by: I.A.C.E.R Srl Via S.Pertini 24/a – 30030**